

**NOROVIRUS GROUP I (RECOMBINANT)
Heat Inactivated Culture Fluid**

PRODUCT DESCRIPTION:

Norovirus Group I (Recombinant) contains ~1Kb sequence of the Norovirus GI genome engineered into a surrogate RNA virus.

Norovirus Group I (Recombinant) is propagated in the Vero cell line.

Heat inactivated culture fluids are sold in 1.0 mL aliquots, and shipped on dry ice.

Each lot of viral culture fluid is assayed for its Tissue Culture Infective Dose (TCID₅₀) before being heat inactivated.

Custom orders are available, including specific volumes and package sizes.

INTENDED USE:

This product is intended for research, product development testing, or quality assurance testing. ***This product can only be used to test for Norovirus in Molecular Assays.***

Heat inactivated culture fluids are sold as consumable testing materials. Applications include:

- Nucleic Acid / Molecular Testing
- Limit of Detection (LOD) Studies
- Cross-reactivity Studies
- Other Viral-based Assays

ETIOLOGIC STATUS/BIOHAZARD TESTING:

Norovirus is a Biosafety Level 2 organism.

Viral inactivation is verified after heat inactivation by the absence of viral growth in validated tissue culture based infectivity assays.

PRECAUTIONS:

USE UNIVERSAL PRECAUTIONS when handling this product!

Heat inactivated culture fluids are treated by a method validated to be effective for virus inactivation. However, no method can be guaranteed 100% effective.

This material should be handled as if capable of transmitting infectious agents.

RECOMMENDED STORAGE:



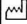
Viral culture fluid should be stored at -65°C or below.

To avoid repeat freeze-thaws, which could negatively impact product performance, store heat inactivated culture fluids in smaller working aliquots.

DO NOT USE IN HUMANS. FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under section 351 of the Public Health Service Act, or for any other product intended for administration to humans.

This product was manufactured in a facility which has a Quality Management System that is ISO 13485 certified.

REF	Catalog Number		Temperature Limitation
LOT	Lot Number		Biological Risk
RUO	For Research Use Only		Date of Manufacture

ZeptoMetrix Corporation
878 Main Street
Buffalo, NY 14202
800-274-5487
www.zeptometrix.com