

**\*HUMAN METAPNEUMOVIRUS (hMPV) 5 Type B1  
Heat Inactivated Culture Fluid**

**PRODUCT DESCRIPTION:**

Human Metapneumovirus is an enveloped virus that contains a negative-sense single-stranded RNA.

hMPV 5 Type B1 (Strain: Peru3-2003) is propagated in the LLC-mk2 cell line.

Heat inactivated culture fluids are sold in 1.0 mL aliquots, and shipped on dry ice.

Each lot of viral culture fluid is assayed for its Tissue Culture Infective Dose (TCID<sub>50</sub>) before being heat inactivated.

Custom orders are available, including specific volumes and package sizes.

**INTENDED USE:**

This product is intended for research, product development testing, or quality assurance testing.

Heat inactivated culture fluids are sold as consumable testing materials. Applications include:

- Nucleic Acid / Molecular Testing
- Cross-reactivity Studies
- Other Viral-based Assays

**ETIOLOGIC STATUS/BIOHAZARD TESTING:**

hMPV is a Biosafety Level 2 organism.

Viral inactivation is verified after heat inactivation by the absence of viral growth in validated tissue culture based infectivity assays.

**PRECAUTIONS:**

USE UNIVERSAL PRECAUTIONS when handling this product!

Heat inactivated culture fluids are treated by a method validated to be effective for virus inactivation. However, no method can be guaranteed 100% effective.

This material should be handled as if capable of transmitting infectious agents.

**RECOMMENDED STORAGE:**

Viral culture fluid should be stored at -65°C or below.





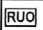

To avoid repeat freeze-thaws, which could negatively impact product performance, store heat inactivated culture fluids in smaller working aliquots.

**DO NOT USE IN HUMANS. FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.**

These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under section 351 of the Public Health Service Act, or for any other product intended for administration to humans.

**\*Note: "This product is sold by ZeptoMetrix Corporation under license from ViroNovative B.V. under patent applications, including U.S. Patent Applications 10/371,099 and 10/371,122, and any patents that issue from applications related to PCT/NL02/00040 and PCT/US03/05271."**

This product was manufactured in a facility which has a Quality Management System that is ISO 13485 certified.

	Catalog Number		Temperature Limitation
	Lot Number		Biological Risk
	For Research Use Only		Date of Manufacture

**ZeptoMetrix Corporation**  
878 Main Street  
Buffalo, NY 14202  
800-274-5487  
www.zeptometrix.com